

Video Article

Protocol for Microplastics Sampling on the Sea Surface and Sample Analysis

Manca Kovač Viršek¹, Andreja Palatinus¹, Špela Koren¹, Monika Peterlin¹, Petra Horvat², Andrej Kržan²

¹Sector for Marine Waters, Institute for water of the Republic of Slovenia

²Laboratory for Polymer Chemistry and Technology, National Institute for Chemistry Slovenia

Correspondence to: Manca Kovač Viršek at manca.virsek@izvrs.si

URL: <http://www.jove.com/video/55161>

DOI: [doi:10.3791/55161](https://doi.org/10.3791/55161)

Keywords: Environmental Sciences, Issue 118, microplastics, protocol, sea surface, microplastics sampling, manta net, microplastics analysis, FTIR spectroscopy, micro FTIR spectroscopy

Date Published: 12/16/2016

Citation: Kovač Viršek, M., Palatinus, A., Koren, Š., Peterlin, M., Horvat, P., Kržan, A. Protocol for Microplastics Sampling on the Sea Surface and Sample Analysis. *J. Vis. Exp.* (118), e55161, doi:10.3791/55161 (2016).

Abstract

Microplastic pollution in the marine environment is a scientific topic that has received increasing attention over the last decade. The majority of scientific publications address microplastic pollution of the sea surface. The protocol below describes the methodology for sampling, sample preparation, separation and chemical identification of microplastic particles. A manta net fixed on an »A frame« attached to the side of the vessel was used for sampling. Microplastic particles caught in the cod end of the net were separated from samples by visual identification and use of stereomicroscopes. Particles were analyzed for their size using an image analysis program and for their chemical structure using ATR-FTIR and micro FTIR spectroscopy. The described protocol is in line with recommendations for microplastics monitoring published by the Marine Strategy Framework Directive (MSFD) Technical Subgroup on Marine Litter. This written protocol with video guide will support the work of researchers that deal with microplastics monitoring all over the world.

Video Link

The video component of this article can be found at <http://www.jove.com/video/55161/>

Introduction

Microplastic pollution in the sea represents a growing concern to contemporary society, due to the constant increase in plastic production and its subsequent disposal and accumulation in the marine environment¹. Even if plastic macro litter would no longer enter the seas, microplastic pollution would continue to grow due to fragmentation of already existing plastic litter in the sea². The majority of microplastic pollution studies were carried out in marine and fresh water ecosystems and mainly addressed sea surface pollution³.

The term microplastic refers to plastic particles smaller than 5 mm in size⁴. This term describes a heterogeneous mixture of particles, which can differ in size (from a few microns to several millimeters), color and shape (from very different shapes of fragments to long fibers). Microplastic particles can be of a primary or secondary origin⁵. Microplastic of primary origin is manufactured as small particles used in the cosmetics industry (pilling crème etc.) or chemical industry as precursor for other plastic products (e.g. plastic pellets used in plastic industry). Microplastic of secondary origin arise via the degradation of larger plastic pieces in the environment due to physical and chemical processes, induced by light, heat, oxygen, water and organisms⁶. In 2015, four types of microplastic sources were defined: larger plastic litter, cleaning products, medicines and textiles⁶. The main source (80 %) of larger plastic litter is assumed to be land based⁷. Microplastic from cosmetic products, medicines and textile enters water ecosystems through sewage and storm waters⁶. Microplastic particles most frequently found in water ecosystems are fragments from larger plastic litter and textile fibers⁸.

Microplastics have several negative effects on the environment. Their small size allows them to enter the food web through ingestion by marine organisms^{9, 10}. Ingested particles can cause physical damage or block the digestive system of animals¹¹. Particles can also be carriers of persistent organic pollutants (POPs). Their hydrophobic surface and favorable ratio of large surface area to small volume, enables POPs to adsorb onto the microplastics¹². In the environment or digestive systems of animals who ingest them, POPs and other plastic additives can be leached from microplastic particles¹³.

Previous studies reported the ubiquitous presence of microplastics in the marine environment³, from the water column to the bottom sediments. The threat of microplastic pollution was already identified by the Marine Strategy Framework Directive in the EU and, consequently, mandatory monitoring of microplastics was advised¹⁴. Accordingly, the EU Technical Subgroup on Marine Litter (TSG-ML) prepared recommendations for monitoring of microplastics in the European seas¹⁵. Thus, the video guidelines for microplastics sampling are of high importance, as they support comparative monitoring and a coherent management process all over the world.

This protocol was developed within the DeFishGear project for the first monitoring of microplastic pollution in the Adriatic Sea. Recommendations from the document "Guidance on Monitoring of Marine Litter in European Seas" by TSG-ML¹⁵ were taken into account. This protocol describes the methodology for microplastics sampling on the sea surface, separation of microplastics from the samples, and chemical analysis of microplastic particles to confirm that particles are from plastic material and to identify the type of plastic. Sampling was done by the use of a

manta net, which is the most suitable equipment for sampling in calm waters¹⁶. Separation of microplastics from the samples was carried out by visual identification using a stereomicroscope. Isolated particles were later chemically identified using Fourier transform infrared (FTIR) spectroscopy and micro FTIR spectroscopy.

Protocol

1. Sampling of microplastics on the sea surface

1. Deploy the manta net from the side of the vessel using a spinnaker boom or »A-frame« using lines and karabiners.
2. Deploy the manta net out of the wake zone (approx. 3 - 4 m distance from the boat) in order to prevent collecting water affected by turbulence inside the wake zone.
3. Write down the initial GPS coordinates and initial time in the data sheet.
4. Start to move in one straight direction with a speed of approx. 2 - 3 knots for 30 min and begin the time measurement.
5. After 30 min stop the boat and write down final GPS coordinates, the length of the route (the most correct way is to calculate the length from the GPS coordinates) and the average boat speed into the data sheet provided and lift the manta net out of the water.
6. Rinse the manta net thoroughly from the outside of the net with seawater using a submersible pump or water from the boat water reservoir. Rinse in the direction from the manta mouth to the cod end in order to concentrate all particles adhered to the net into the cod end.
Note: Never rinse the sample through the opening of the net in order to prevent contamination.
7. Safely remove the cod end and sieve the sample in the cod end through a 300 µm mesh size sieve or less.
8. Rinse the cod end thoroughly from the outside and pour the rest of the sample through the sieve. Repeat this step until there are no longer any particles inside the cod end.
9. Concentrate all material on the sieve in one part of the sieve.
10. With the use of a funnel, rinse the sieve into a glass jar or plastic bottle using 70 % ethanol.
11. Close the bottle, wipe it with paper towels and label the lid and outside of the jar with the sample name and date with waterproof marker (you should also put a second label written with a pencil on velum paper in a jar to avoid the possible loss of the sample name due to the erased label on the jar). Transfer labeled plastic bottle into the cool box.
Note to general sampling conditions: The wind speed should not be more than 2 Beaufort, since the waves are too high and the net is not stable on the sea surface. It is important to maintain a steady linear course at a constant speed during the trawls. Half of the manta net opening should be submersed during sampling. Duration of sampling should be 30 min (in cases where there is a large amount of natural material, e.g. plankton bloom, the duration of sampling can be shorter). Avoid the use of plastic tools and containers. Avoid synthetic clothing (e.g. fleece), ropes and contact of manta net with vessel to prevent contamination of the sample. Be very careful not to damage the manta net or the boat hull while deploying and capturing the net.

2. Separation of microplastics from the sea surface samples

1. If the sample does not contain any items larger than 25 mm and appears to be clean, continue directly with step 3.
2. Pour sample through the sieve (≤300 µm mesh size) and remove all natural or artificial litter objects of a size >5 mm (macro and mezzo litter) from the sample, using visual identification and tweezers. Be careful to rinse each removed object carefully with distilled water in order to remove any microplastic litter adhered to it. Store all natural and artificial litter objects in separate containers. Dry all natural and artificial litter objects in a desiccator (or in the open air, but in a closed dish) and weigh them. Identify all litter objects >25 mm (macro litter) according to the Master List of Categories of Litter Items¹⁶.
3. After removing all larger objects, concentrate all remaining pieces in one part of the sieve using squirt bottles or tap water. Pour the sample into a glass container using a minimum amount of 70 % ethanol with the help of a funnel.
Note: In this step the use of 70 % ethanol is crucial to preserve the sample. Also in the step of visual inspection of the sample, ethanol helps to discolor the organisms and colorful plastics therefore become easier to find.
4. Take a small amount of the sample (subsample) and pour it into a glass Petri dish. Analyze the sample with the use of a stereomicroscope (20 - 80x zoom) and search for microplastic particles.
5. Each microplastic particle should be categorized into one of the categories listed in Table 1 and put into a Petri dish or other glass vials, marked with a category name. The Petri dish needs to be closed at all times.
Note: When separating microplastics from your sample be conservative and select more rather than less particles for the analysis. The real chemical structure of particles will still be determined later on. Be sure to analyze larger objects from all sides as microplastics may be stuck and therefore hidden under larger items. It may also be helpful to move already analyzed objects to one side of the Petri dish.
6. Put the Petri dish under the microscope with measuring equipment (ocular ruler calibrated by the micrometer slide or image analysis software) and measure the size of each particle (measure the longest diagonal), except filaments, and note its color. Each subsample should be reviewed by another person. Be careful to rinse the glass container containing the sample so that all particles adhering to the glass walls are washed into the Petri dish.
7. Weigh the microplastic particles of each category separately by the use of analytical scale. Microplastic particles need to be dried prior to weighing. The closed Petri dish can be put in a desiccator or the samples can be left to dry in a closed dish till particles became dry (the weight of closed petri dish with particles is constant).
8. Identify micro litter. When analyzing a sample in search of microplastics, please consider that some particles will be easily visible (color, shape, size) while others may be trickier to find. Below are a few features that identify microplastic particles in the sample: For example, no cell structure, uneven, sharp, crooked edges, uniform thickness, distinctive colors (blue, green, yellow, etc.).

3. Chemical identification of microplastics

1. **ATR-FTIR spectroscopy**
 1. Prior to the analysis clean the detection system with alcohol and a lint free cloth.

2. Record a background spectrum. Place the sample on the sample holder and collect the spectra. Identify the obtained ATR- FTIR spectra using an automated comparison of the obtained spectrum with spectra in a database.

2. Micro ATR-FTIR spectroscopy

1. Prior to the analysis clean the detection system with alcohol and a lint free cloth.
2. Place the sample on a glass filter. Note: Other filters can be used but their polymer nature can interfere with the characterization.
3. Place the filter with the sample on the automatic scanning table and use the joystick to locate the sample.
4. Record an optical image and mark an area (e.g. 20 by 20 μm) where the sample will be characterized.
5. Record a background spectrum.
6. Place the sample on the sample holder and collect the spectra at the predefined location.
7. Identify the obtained micro ATR-FTIR spectra using an automated comparison of the obtained spectrum with spectra in a database.

Representative Results

The first result of the described protocol are microplastic particles categorized into six categories according to their visual features (Table 1). The first category, and usually the most abundant one, are fragments (Figure 1). They are rigid, thick, with sharp crooked edges and an irregular shape. They can be in a variety of different colors. The second category are films (Figure 2). They also appear in irregular shapes, but in comparison with fragments, they are thin and flexible and usually transparent. The third category are pellets (Figure 3), usually originating from the plastics industry. They are irregular, round shapes, and normally bigger in size, around 5 mm in diameter. They are usually flat on one side and can be of various colors. The fourth category are granules (Figure 4). In comparison with pellets, they have a regular round shape and usually a smaller size, around 1 mm in diameter. They appear in natural colors (white, beige, brown). The fifth category are filaments (Figure 5). They are, next to fragments, the most abundant type of microplastic particles. They can be short or long, with different thicknesses and colors. The last category are foams (Figure 6). They most often come from large particles of styrofoam. They are a soft, irregular shape and white to yellow in color.

The main result of microplastics sampling and sample analysis is the number of microplastic particles per sample. These data can be further normalized per km^2 . The formula used for normalization is:

microplastic particles per sample / sampling area,

where sampling area is calculated by multiplying sampling distance by the width of the opening of the manta net (Tables 2, 3; Figure 7). In addition, particles can be analyzed with image analysis software. The results include maximum length and area of each particle (Table 4). Figure 8a show particles before image analysis and Figure 8b is after image analysis, where each particle is measured and numbered. Lastly, a chemical analysis of the total or highest possible number of particles per sample is recommended. Using Fourier transform infrared spectroscopy a spectrum of the selected particle is acquired, as shown on Figure 9. This spectrum is then compared with the spectra from the software library (Figure 10). The final result will show if a given particle is plastic or not and indicate the type of plastic from the chemical structure.

1	Fragments
2	Films
3	Pellets
4	Granules
5	Filaments
6	Foams

Table 1: Categories of microplastic particles.



Figure 1: Example of particles from category: Fragments. [Please click here to view a larger version of this figure.](#)

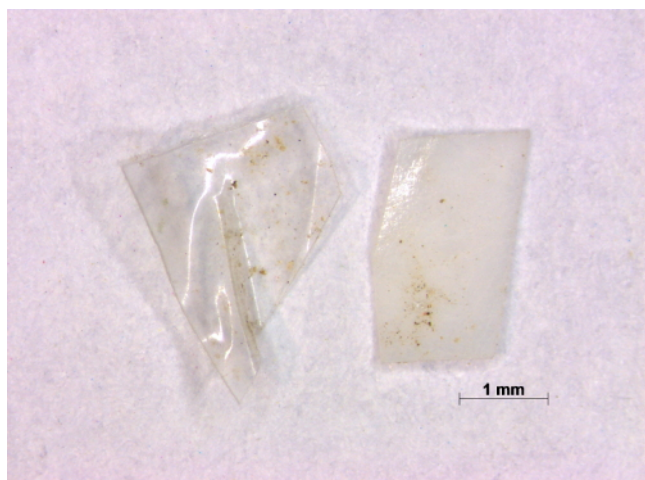


Figure 2: Example of particles from category: Films. [Please click here to view a larger version of this figure.](#)



Figure 3: Example of particles from category: Pellets. [Please click here to view a larger version of this figure.](#)

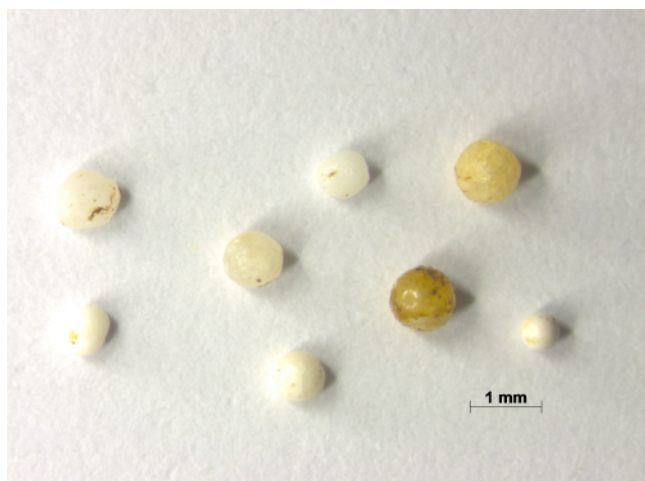


Figure 4: Example of particles from category: Granules. [Please click here to view a larger version of this figure.](#)

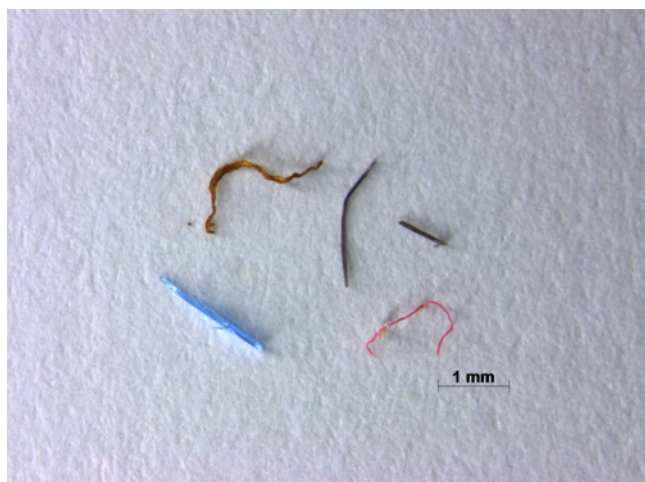


Figure 5: Example of particles from category: Filaments. [Please click here to view a larger version of this figure.](#)

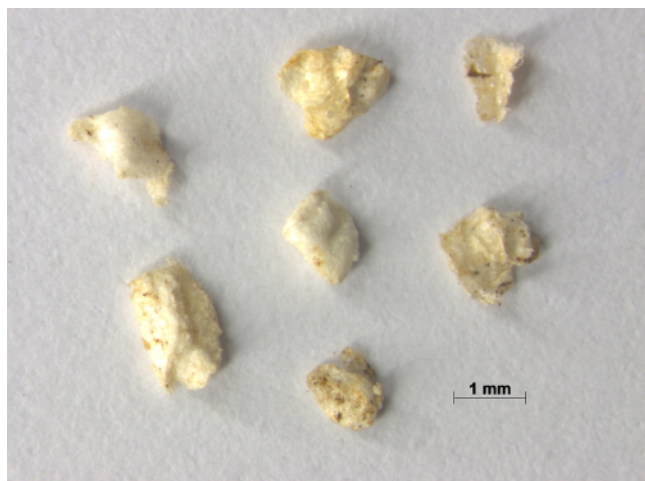


Figure 6: Example of particles from category: Foams. [Please click here to view a larger version of this figure.](#)

Sampling distance [km]	2
Manta width [km]	0.0006
Sampling area [km ²]	0.0012

Table 2: Example of data from survey, used for calculation of microplastic particles per km².

	No	No/km ²
fragments	301	250833
films	45	37500
pellets	15	12500
granules	8	6667
foams	33	27500
filaments	223	185833

Table 3: Example of results from survey, where the categorized data into 6 groups are counted and normalized per km² (No - number of particles).

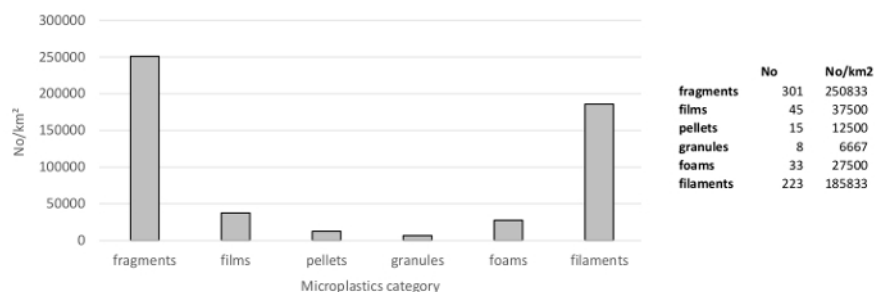


Figure 7: Example of representative results after visual categorization of particles (No - number of particles). [Please click here to view a larger version of this figure.](#)

Index Region	Area [mm ²]	Maximum length [mm]
1	8.010	5.506
2	10.517	5.628
3	12.185	5.429
4	3.367	3.367
5	2.475	2.155
6	1.809	2.943
7	6.604	5.238
8	5.779	4.037
9	4.472	3.791
10	16.907	5.355
11	7.246	3.733
12	7.867	4.622
13	6.411	5.056
14	3.281	3.070
15	12.937	5.554
16	6.709	3.716

Table 4: Example of image analysis results where area [mm²] and maximum length [mm] of each particle are measured.

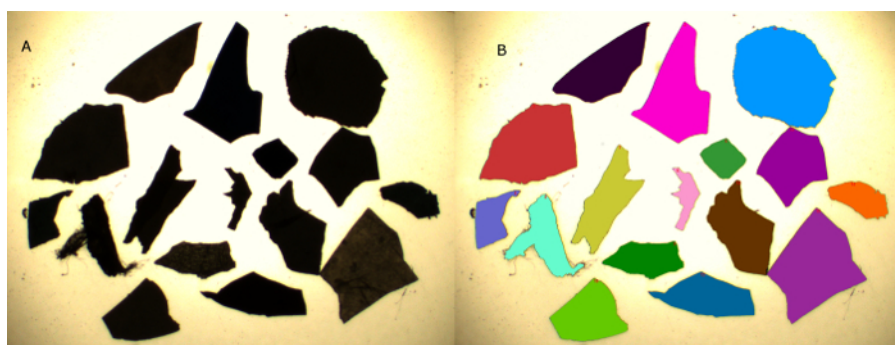


Figure 8: Example of image acquired a) before and b) after image analysis of particles with image analysis software. [Please click here to view a larger version of this figure.](#)

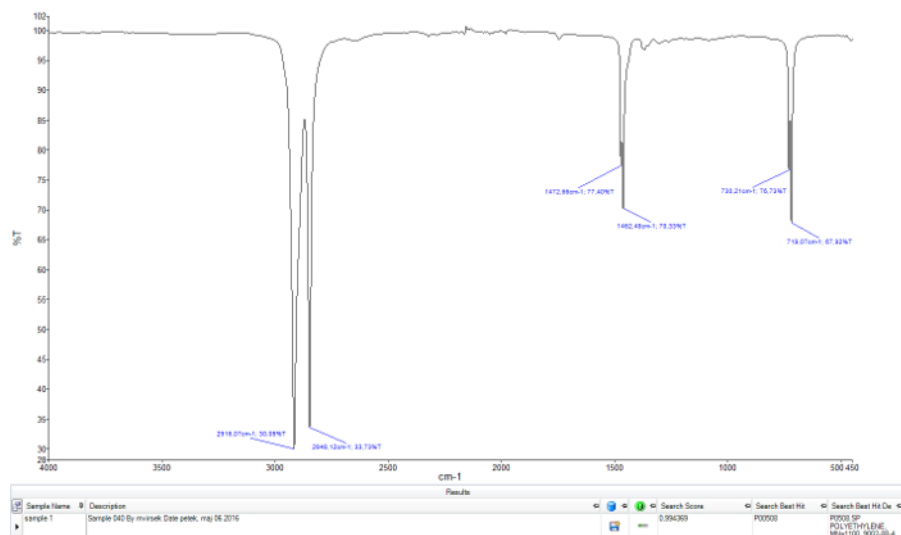


Figure 9: Example of a spectra measured on a selected particle with marked peaks and their wavelengths [cm^{-1}]. Please click here to view a larger version of this figure.

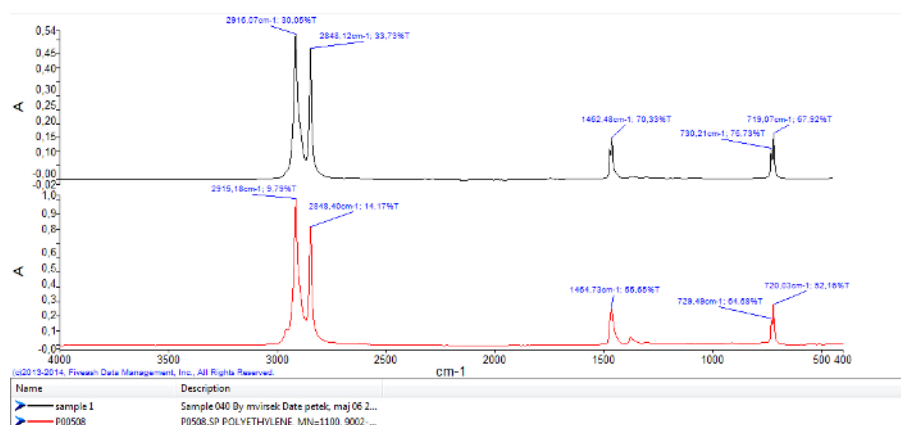


Figure 10: Example of comparison of acquired spectra from selected particle to best match from the ATR-FTIR spectra library. Please click here to view a larger version of this figure.

Discussion

Microplastics sampling on the sea surface by manta net is a widely used method for the sampling of microplastics on the sea surface, but to date there has been no unified methodology. A large volume of water can be filtered through the manta net, thus the possibility of trapping a relevant number of microplastics is high and the results are perceived to be reliable. Comparability of results among different samples is assured by normalization. In our case, the concentrations were related to the sampled area by multiplying trawl distance by the horizontal width of the net opening. Another option is to use a flow meter, fixed at the net opening. The use of a flow meter is possible since the manta net with its lateral wings is very stable on the sea surface and therefore hopping on the waves is minimal. A flow meter records the volume of filtered water and thus enables the normalization of results per volume of sampled water¹⁶.

The most frequently used manta nets have around 300 μm mesh size and are 3 - 4.5 m long. These dimensions were optimized to avoid clogging of the net and to allow the sampling a volume of water as large as possible. Trawling speed is recommended to be between 2 - 3 knots, but it is dependent on wave height, wind speed and sea currents. It is very important that the manta net is under supervision the whole time during sampling and if it starts hopping, the trawling speed must be reduced. The trawling time is recommended to be around 30 min, but depends on seston concentrations. It can happen that seston sometimes clogs the manta net. In this case the trawling has to be stopped immediately, otherwise the microplastic particles can be lost and the net can get damaged. Manta net is the most often fixed from the side of the vessel. This is also the most suitable option, while the manta net is surely out of the wake zone. In some surveys manta net was fixed from the stern of the vessel^{17, 18}, but in that case you have to be sure that the net is out of the wake zone. The distance, on which the trawl is set for sampling, should be determined individually, since the zone of turbulences caused by the vessel varies from the size of the vessel and from the speed of the boat^{19, 20}.

Separation of microplastic particles from the sea surface samples is most often done just by visual identification²¹. Particles bigger than 1 mm can be identified easily by the naked eye, while particles smaller than 1 mm require the use of a stereomicroscope. To reduce the possibility of confusing the non-plastic particles with plastic ones, using the polarization light on stereomicroscopes is recommended. The possibility of

misidentification of plastic particles gets higher with smaller particles. Thus particles >0.5 mm can only be identified visually²¹, by the use of stereomicroscope. For particles smaller than 0.5 mm an additional, more accurate method is required e.g. micro ATR-FTIR spectroscopy²¹.

During the process of microplastics separation from the sample the possibility of sample contamination with the airborne filaments is very high. For this reason, control Petri dishes left open on the working table are strongly recommended for the identification of potential contaminant airborne particles. Namely, the quality of the data strongly depends on: 1) the precision of the person working with the sample, 2) the quality and magnification of the stereomicroscope, and 3) the quantity of organic matter in the sample¹⁶. After visual identification it is strongly recommended to analyze the sorted particles with one of the available techniques for chemical identification of the material⁸.

Several methods exist for polymer identification, among which the FTIR spectroscopy and Raman spectroscopy are the most frequently used²². FTIR and Raman spectroscopy are complementary techniques and their accuracy is similar. In our protocol, the FTIR and micro FTIR spectroscopy with "attenuated total reflectance" (ATR) are presented. They are simple to use and they enable fast and accurate results. Plastic polymers possess highly specific infrared (IR) spectra with distinct band patterns, thus making IR spectroscopy an optimal technique for the identification of microplastics²¹. The energy of IR radiation excites a specific molecular vibration when interacting with a sample, which enables the measurement of characteristic IR spectra²². FTIR spectroscopy can also provide additional information on particles, such as intensity of oxidation²³ and level of degradation²⁴. While ATR-FTIR is suitable for chemical identification of larger particles (>0.5 mm), micro ATR-FTIR spectroscopy can provide information on the chemical structure of particles <0.5 mm, as it combines the function of a microscope and an infrared spectrometer.

Before using FTIR and micro FTIR spectroscopy, microplastic particles have to be previously dried, since water strongly absorbs IR radiation²², and purified, in case they are covered with biofilms and/or other organic and inorganic adherents, which can influence the IR spectra. The most non-invasive way to purify samples is by stirring and rinsing with fresh water²⁵. If this is not enough, then the use of 30 % hydrogen peroxide is recommended. All other methods can have negative effects on the microplastic particles (e.g. ultrasonic cleaning can further break particles, strong acidic or alkaline solutions can damage several plastic polymers, etc.) and therefore their use is not recommended. More promising is the use of a sequential enzymatic digestion as a plastic friendly purification step. Purification using different technical enzymes (e.g. lipase, amylase, proteinase, chitinase, cellulase, proteinase-K) has been successfully applied to reducing a biological matrix of plankton and thus proved to be a valuable technique to minimize matrix artifacts during FTIR spectroscopy measurements²².

Separation of microplastics by visual identification and chemical identification of selected particles are both extremely time-consuming processes. This work has to be done by an accurate and patient person who has experience with stereomicroscopes, not only in recognizing the plastic particles, but also in recognizing biological matter. Even an experienced person cannot discriminate all potential microplastic particles unambiguously from chitin or diatom fragments²². Therefore, the error rate of visual sorting ranges from 20 %²⁶ to 70 %²¹ and increases with decreasing particle size.

Disclosures

The authors have nothing to disclose.

Acknowledgements

The development of this protocol was founded by IPA Adriatic Cross-border Cooperation Programme 2007-2013, within the DeFishGear project (1°str/00010).

References

1. Law, K. L., *et al.* Plastic accumulation in the North Atlantic subtropical gyre. *Science*. 329 (5996), 1185-1188 (2010).
2. Thompson, R. C. Microplastics in the marine environment: Sources, consequences and solutions. In *Marine anthropogenic litter*. Springer International Publishing. 185-200 (2015).
3. Lusher, A. Microplastics in the marine environment: distribution, interactions and effects. In: *Marine anthropogenic litter*. Springer International Publishing. 245-307 (2015).
4. Arthur, C., Baker, J., & Bamford, H. Proceedings of the International Research Workshop on the Occurrence, Effects, and Fate of Microplastic Marine Debris, September 9-11, NOAA Technical Memorandum NOS-OR&R30. (2008).
5. Andrady, A.L. Microplastics in the marine environment. *Marine pollution bulletin*. 62(8), 1596-1605 (2011).
6. Browne, M.A. Sources and pathways of microplastics to habitats. In: *Marine anthropogenic litter*. Springer International Publishing. 229-244 (2015).
7. UNEP's REGIONAL SEAS PROGRAMME. *Marine litter: an analytical overview*. UNEP, (2005).
8. Wal, M., *et al.* SFRA0025: Identification and Assessment of Riverine Input of (Marine) Litter, *Final Report for the European Commission DG Environment under Framework Contract No ENV.D.2/FRA/2012/0025*. (2015).
9. Setälä, O., Fleming-Lehtinen, V., Lehtiniemi, M. Ingestion and transfer of microplastics in the planktonic food web. *Environmental pollution*. 185: 77-83 (2014).
10. Farrell, P., Nelson, K. Trophic level transfer of microplastic: *Mytilus edulis*. (L.) to *Carcinus maenas*. (L.). *Environmental Pollution*. 177, 1-3 (2013).
11. Wright, S.L., Thompson, R.C., Galloway, T.S. The physical impacts of microplastics on marine organisms: a review. *Environmental Pollution*. 178: 483-492 (2013).
12. Bakir, A., Rowland, S.J., Thompson, R.C. Transport of persistent organic pollutants by microplastics in estuarine conditions. *Estuarine, Coastal and Shelf Science*. 140, 14-21 (2014).
13. Cole, M., Lindeque, P., Halsband, C., Galloway, T.S. Microplastics as contaminants in the marine environment: a review. *Marine pollution bulletin*. 62(12), 2588-2597 (2011).

14. Zarfl, C., et al. Microplastics in oceans. *Marine Pollution Bulletin*. 62, 1589-1591 (2011).
15. Hanke, G., et al. MSFD GES technical subgroup on marine litter. Guidance on monitoring of marine litter in European Seas. *Luxembourg: Joint Research Centre-Institute for Environment and Sustainability*. Publications Office of the European Union. (2013).
16. Gerdts, G. Methodology used for the detection and identification of microplastics - A critical appraisal. In: *Marine anthropogenic litter*. Springer International Publishing. 201 - 227 (2015).
17. Kang, J. H., Kwon, O. Y., Lee, K. W., Song, Y. K., Shim, W. J. Marine neustonic microplastics around the southeastern coast of Korea. *Marine pollution bulletin*. 96(1), 304-312 (2015).
18. Lusher, A. L., Tirelli, V., O'Connor, I., Officer, R. Microplastics in Arctic polar waters: the first reported values of particles in surface and sub-surface samples. *Scientific reports*. 5 (2015).
19. Shu, J.-J. Transient Marangoni waves due to impulsive motion of a submerged body. *International Applied Mechanics*. 40 (6): 709-714 (2004).
20. Rabaud, M., Moisy, F. Ship wakes: Kelvin or Mach angle?, *Physical Review Letters*. 110 (21): 214503 (2013).
21. Hidalgo-Ruz, V., Gutow, L., Thompson, R. C., Thiel, M. Microplastics in the marine environment: a review of the methods used for identification and quantification. *Environmental science & technology*. 46(6), 3060-3075 (2012).
22. Kuczera, M., Mintenig, S., Lorenz, C., Gerdts, G. Focal plane array detector-based micro-Fourier-transform infrared imaging for the analysis of microplastics in environmental samples. *Environmental Chemistry*. 12(5), 563-581 P. L., Biesinger, M. C., Grifi, M. Plastics and beaches: A degrading relationship. *Marine Pollution Bulletin*. 58(1), 80-84 (2009).
23. Ioakeimidis, C., et al. The degradation potential of PET bottles in the marine environment: An ATR-FTIR based approach. *Scientific reports*. 6, 23501 (2016).
24. McDermid, K.J., McMullen, T.L. Quantitative analysis of small-plastic debris on beaches in the Hawaiian archipelago. *Marine pollution bulletin*. 48, (7), 790-794 (2004).
25. Eriksen, M., et al. Microplastic pollution in the surface waters of the Laurentian Great Lakes. *Marine pollution bulletin*. 77(1-2), 177-182 (2013).