

Figure S5

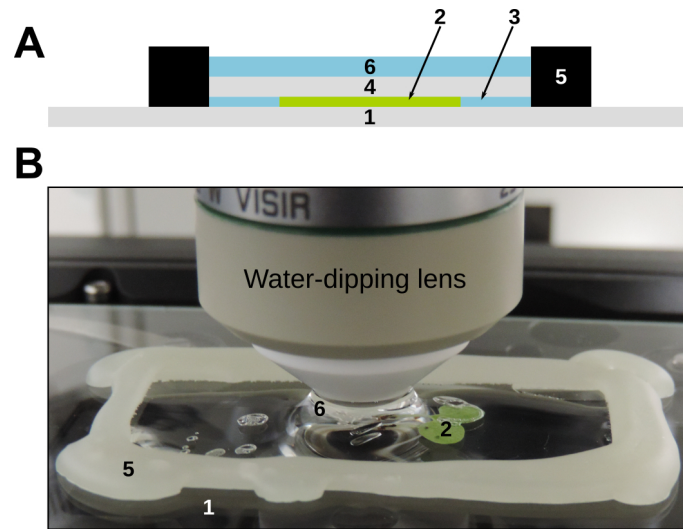


Figure S5. Experimental set-up used for time-lapse imaging of infected *N. benthamiana* roots. **(A)** Schematic representation of the experimental set-up. A one-week-old *N. benthamiana* seedling (2) is mounted in a liquid compartment (3) between a slide (1) and a coverslip (4). The edges of the coverslip are sealed with a 4:1 mix of paraffin and lanolin (5). The upper part of the coverslip is filled with water for use with a water-dipping objective. **(B)** Representative image of the experimental set-up. Numbers refer to the same elements as described previously.