Supplementary Information

Scrodentoids H&I, a pair of natural epimerides from *Scrophularia dentata*, inhibits inflammation through JNK-STAT3 axis in THP-1 cells

Gaohui Mao^{1,2}, Liqin Sun^{1,2}, Jinwen Xu¹, Yiming Li², Ciren Dunzhu³, Liuqiang Zhang²*, Fei Qian¹*

Email: 04100217@163.com (L.Z.); qianfei0517@126.com (F. Q.)

This PDF file includes:

Supplementary text Figures. S1, S2

SI Materials and Methods

Materials

p-JAK2 antibody waspurchased from Abcam, JAk2 antibody was purchased from Cell Signaling Technology(USA). CST PC was a cell extract which is a positive control of p - JAK2 from CST Company.

Method

THP-1 cells were pretreated with SHI for 1 hr, and then stimulated with for 30 min. After the treatments, cell lysates were prepared using radio-immunoprecipitation assay (RIPA) (Beyotime, China) lysis buffer containing protease and phosphatase inhibitor cocktails. The protein concentration was measured using the bicinchoninic acid (BCA) protein assay kit (Beyotime, China). Equal concentrations of protein were separated using a 10% SDS-PAGE and transferred to PVDF membranes. After blocking, the membranes were incubated with antibodies against p-JAK2, JAK2, GAPDH. The membranes were incubated with horseradish peroxidase (HRP-conjugated secondary antibodies (Cell Signaling Technology, USA) and developed using the enhanced chemiluminescence (ECL) detection system (Millpore, USA).Protein Eschborn, Germany) activity was determined using Chemiluminescent Im aging System (Tanon 5200 Multi) and Gel Image System (Tanon, China).

SI Figure legends

Figure. S1 The effect of SHI and SP600125 on p-JAK2/JAK2. Cells were pre-treated with various concentrations of SHI or SP600125 for 1 h following treatment of LPS (1 μ g/mL) for 30 min. Figure A, B, C are presented three independent experiments.

Figure. S2 The effect of CST PC and our samples on p-JAK2. CST PC was a positive control cell extracts from CST Company. Another four samples, THP-1 cells were pre-treated without or with SP600125 for 1 h following treatment of LPS (1

 $\mu g/mL$) for 30 min.

Figures

Figure S1

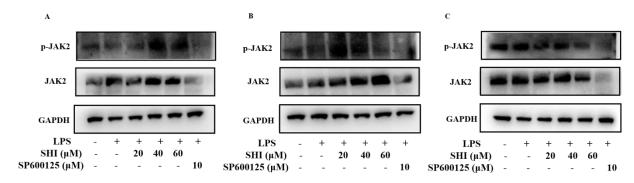


Figure S2

