

## **BBL** **MEDIA for FUNGI**

**In Dehydrated Form or  
Sterilized and Ready for Use**

**LITTMAN MEDIUM  
MYCOPHIL AGAR  
MYCOPHIL BROTH  
SABOURAUD DEXTROSE AGAR  
SABOURAUD MALTOSSE AGAR  
CORN MEAL AGAR  
without added Dextrose**

*Write for Folder No. 3095*

**BALTIMORE  
BIOLOGICAL LABORATORY**

1640 Gorsuch Avenue      Baltimore-18, Md.

## **Reconditioned Microscopes** **ALL STANDARD MAKES**

**B & L — Spencer — Leitz — Zeiss**

### *Specials*

Spencer Vert. Binc. Research  
B & L JK Wide Field  
B & L SAK Wide Field  
B & L Polarizing  
Spencer Model 860 Microtome  
B & L Freezing Clinical  
Spencer & BL Rotary &  
Student Microtomes

*Several Wide field Stereoscopic Type*

*All instruments Guaranteed*

**Complete Repair Service**

**MONROE MICROSCOPE SERVICE**


P. O. Box 656  
Rochester 2, N.Y.

## **LKB PAPER ELECTROPHORESIS EQUIPMENT**

Paper electrophoresis has proved to be a valuable aid in clinical laboratory work. Analyses of serum protein are obtained rapidly and with an accuracy which is sufficient for diagnostic purposes.

The LKB Paper Electrophoresis Apparatus is easy to operate, affords full protection against electric shocks, and has the exceptionally large capacity of 12 simultaneous analyses. The Power Supply is equipped with a timer which disconnects the current after a predetermined running time. The maximum output is 50 ma. at 275 volts. Adaptable to 110, 127, 220, and 240 volts A.C., 50-60 cycles.

*Write for free bulletin L-69*

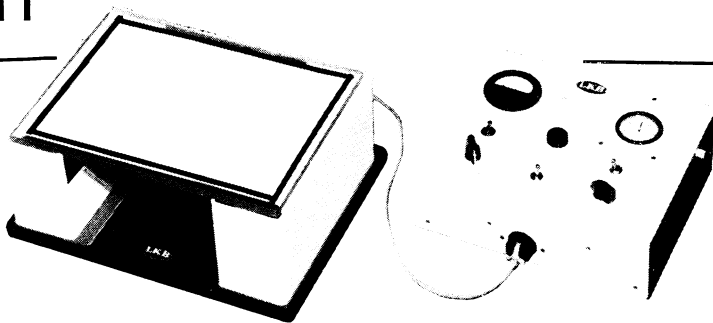

# **Ivan Sorvall, inc.**

210 Fifth Ave.  
New York 10, N.Y.


(In writing to advertisers, please mention the journal—it helps.)

*Improved*  
Model

**MODEL PR-1**

# International Refrigerated Centrifuge

The New International Model PR-1 Portable Refrigerated Centrifuge embodies the many time-proven features found in previous models and, in addition, incorporates important engineering improvements. The result is superior performance in the laboratory without increase in cost.

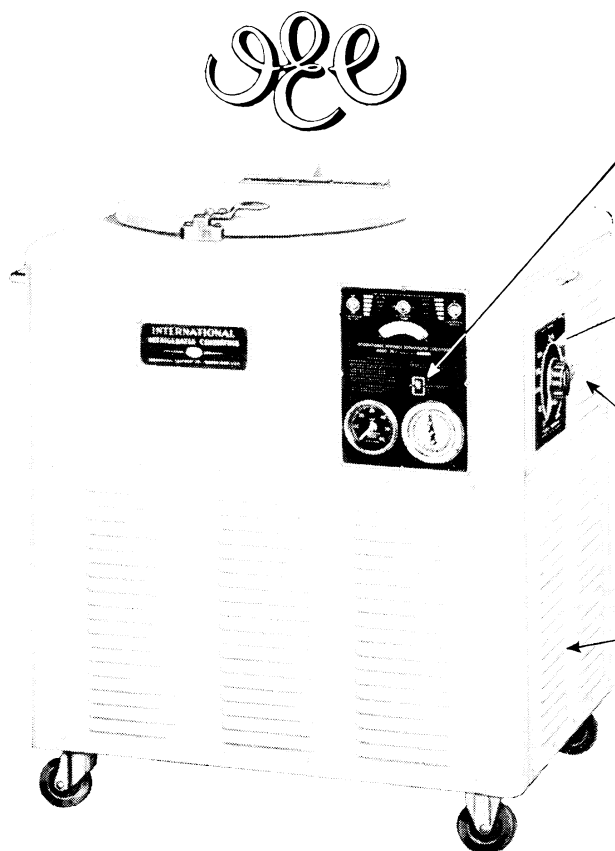

A new type of temperature controller enables the user to pre-select the operating temperature within the Centrifuge chamber. Temperature is controlled within  $\pm 1^{\circ}\text{C}$ . and indicated on a calibrated dial.

An autotransformer speed controller is provided with alternating current models, eliminating heating associated with the rheostat type of controller, and providing uniform, stepless speed control throughout the entire range.

A new hinged panel provides easy access to all instruments and electric controls, thus facilitating maintenance.

A hermetically sealed unit is substituted for the open type compressor heretofore used. This eliminates trouble some vibration caused by the action of the compressor and minimizes the possibility of gas leakage.

Interchangeable accessories are available in wider variety, enhancing the usefulness of the Centrifuge.

Send today for new descriptive Bulletin RC-1952 containing complete details and current prices.

**INTERNATIONAL EQUIPMENT COMPANY**

1284 SOLDIERS FIELD ROAD, BOSTON 35, MASS.



## **BACTO-TRYPSIN 1% for Hemagglutination** **BACTO-HEMAGGLUTINATION BUFFER**

THE USE of trypsinized red cells in Rh determinations has improved both the sensitivity and specificity of the test making the direct procedure as reliable as the more involved indirect antiglobulin tests. Trypsin treated erythrocytes favorably influence the determination of the A, B, S, Le<sup>a</sup>, Le<sup>b</sup> and P factors as well as provide a more sensitive cell suspension for virus hemagglutination studies.

### ***Bacto-Trypsin 1% for Hemagglutination***

is a carefully standardized soluble Trypsin for conditioning erythrocytes for Rh determinations and other hemagglutination tests. Bacto-Trypsin 1% for Hemagglutination is reconstituted by the addition of Bacto-Hemagglutination Buffer solution, yielding a sparklingly clear solution equivalent in tryptic activity to a 1 per cent solution of Trypsin, Difco 1:250.

### ***Bacto-Hemagglutination Buffer***

rehydrated with distilled water is recommended for the preparation of isotonic buffered saline for use in hemagglutination studies, especially for reconstituting Bacto-Trypsin 1% for Hemagglutination and for suspending and washing erythrocytes preparatory for hemagglutination.

***These reagents*** are supplied in companion packages containing 6 vials (10 ml.) of Bacto-Trypsin 1% and 6 vials of Bacto-Hemagglutination Buffer each sufficient for preparing 1000 ml. of buffered saline.

THE DIFCO MANUAL, NINTH EDITION,  
*including descriptions of these media and their use,*  
*is available on request.*

**DIFCO LABORATORIES**  
**DETROIT 1, MICHIGAN**