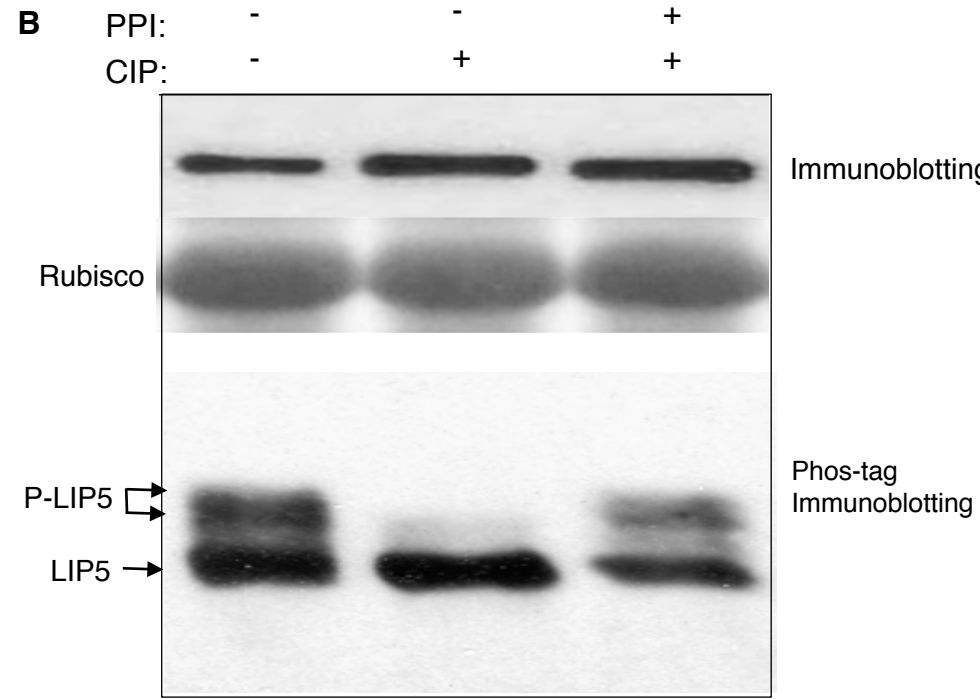
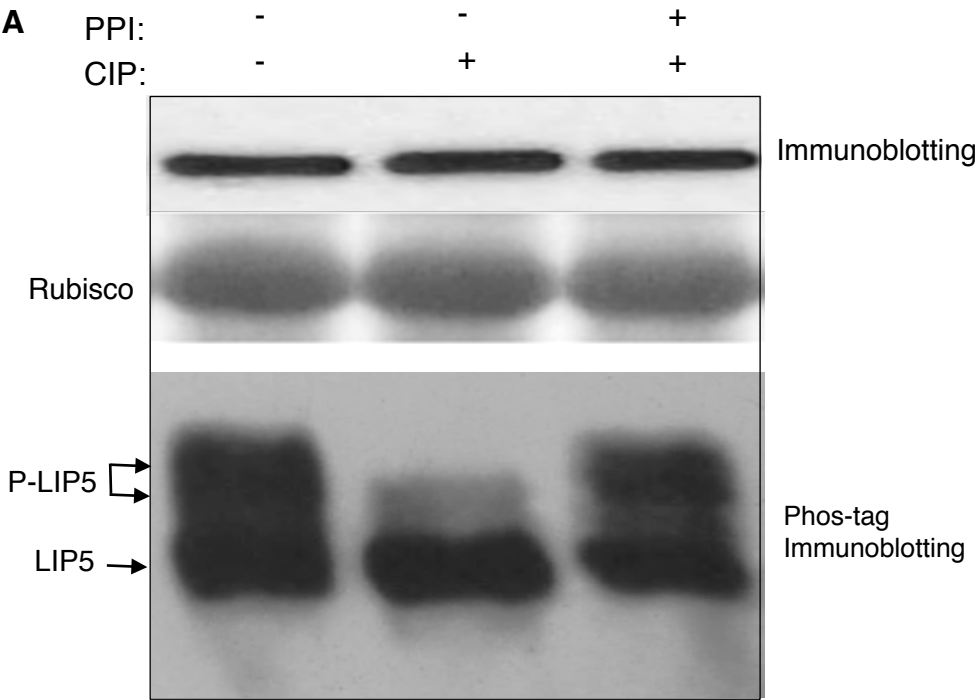


Figure S7



**Figure S7.** Dephosphorylation of *in vivo* Phosphorylated LIP5 Proteins.

Protein extracts were isolated from transgenic *NtMEK2<sup>DD</sup>/myc-LIP5<sup>WT</sup>* at 24 hours after DEX treatment (**A**) or *lip5-1/myc-LIP5<sup>WT</sup>* (**B**) plants at 24 hpi of *Pst*DC3000. The protein extracts were treated at 37°C for 45 minutes with calf intestinal alkaline phosphatase (CIP) in the absence or presence of a phosphatase inhibitor cocktail (10 mM NaF, 7 mM  $\beta$ -glycerophosphate and 5 mM Na-pyrophosphate). Reactions without CIP and phosphatase inhibitors (-) were used as control. The protein extracts were subsequently separated on the regular SDS-PAGE and Phos-tag gels for immunoblot analysis using an anti-myc monoclonal antibody. Rubisco staining of the regular SDS-PAGE gel was used for assessing equal protein loading.