


New ELECTRONICALLY Controlled Water Bath Shaker!

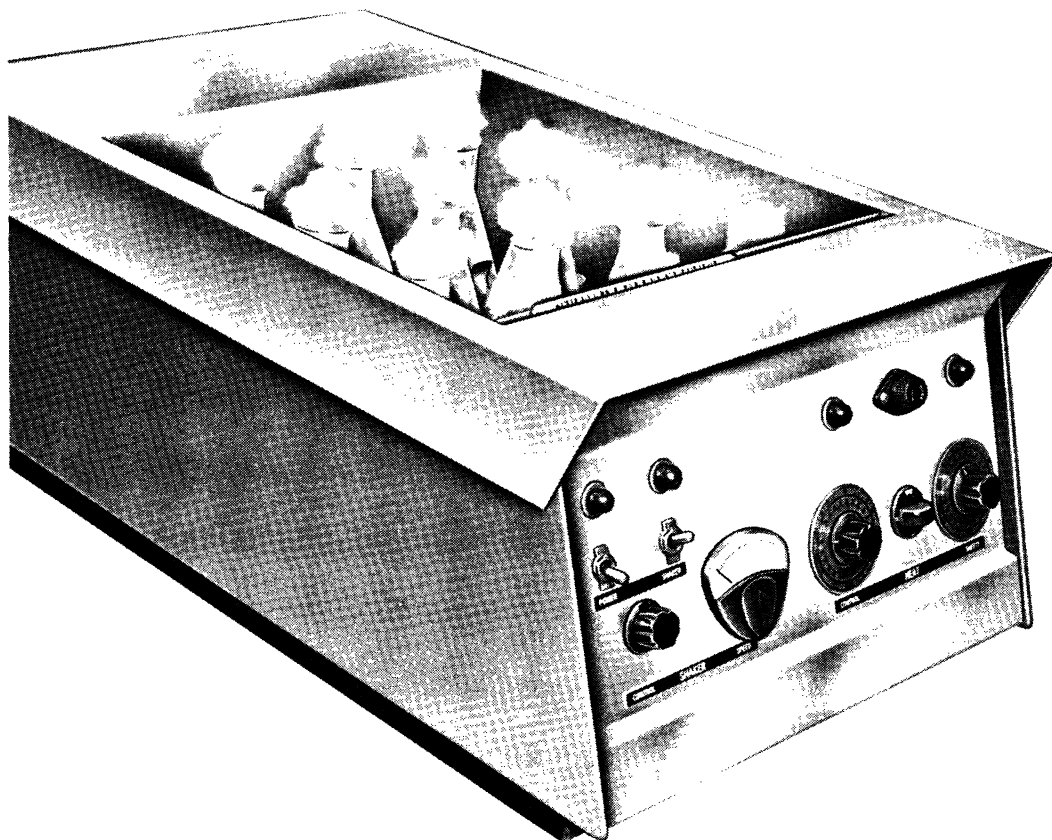

Controls Temperature ELECTRONICALLY within $\pm 0.25^{\circ} \text{C}$!

Controls Agitation ELECTRONICALLY from 40-400 rpm!

Controls Water Level ELECTRONICALLY!

Controls Heating and Tap-Water Cooling ELECTRONICALLY!

Other new engineering features: Auxiliary safety thermostat protects against overheating. Safety lock-knobs prevent accidental change of thermostat settings. Right-angle mercury thermometer is designed for easier reading without removal, and is also protected against breakage. Extra large shaker capacity.

Send for catalog G 77B/65



New Brunswick Scientific Company, Inc.
1130 Somerset Street, New Brunswick, N. J. 08903

12 advantages with new FALCON sterile disposable Plastic Pipets

sterile / disposable / inert / safe / cotton-plugged / chemically
clean / convenient / accurate / space-saving /
non-wettable / uniform / easy-to-read

new FALCONTM rigid plastic disposable pipet

Sterile, plugged, packaged, identified—always ready to use—these new pipets reduce the ever-increasing laboratory work load. You will recognize the speed and convenience possible in all pipetting procedures where uniformity and sterility are mandatory. The non-wettable feature significantly reduces drainage "hold-back" errors and provides an added bonus—the meniscus is minimized, making easy, rapid readings a routine. Unusual clarity of markings is further assured for every test by the use of a fresh pipet. Available: 1.0 and 1.1 ml. pipets in plastic packs of 25; also 1.0 ml. pipets (0.01 or 0.1 ml. div.) individually wrapped, in boxes of 100.

Products of B-D LABORATORIES are available through your local distributor. Further information on request.



B-D LABORATORIES INC., RUTHERFORD, NEW JERSEY



Becton
Dickinson
Laboratory



Falcon
Plastics



Science
Education
Products

In Canada: Becton, Dickinson & Co., Canada, Ltd., Clarkson, Ontario
Overseas: Becton, Dickinson & Co., S.A., P.O. Box 1173, Colon Free Zone, Panama

FALCON, B-B-L, and B-B-L ARE TRADEMARKS © 1971



Antimicrobial Agents and Chemotherapy 1964

Proceedings of the
Fourth Interscience
Conference on
Antimicrobial Agents
and Chemotherapy

NEW YORK, N.Y.

October 26–28, 1964

Sponsored by the American Society for Microbiology

Reports of basic research on the new antibiotics almarcetin, enteromycin carboxamide, everninomicin, hamycin, LL-AL471E, LL-AP191, monicamycin, rubiflavin, U-11,973, U-21,799, U-13,714, U-11,092, U-22,956, and others

Clinical evaluation of caphalothin, nafcillin, ampicillin, oxacillin, cloxacillin, dicloxacillin, methicillin, lincomycin, spectinomycin, nalidixic acid, tolnaftate, nitrofurantoin, colistin, capreomycin, gentamicin, rifomycin, and others

Symposia on Penicillins and Cephalosporins, Synthetic Antimicrobials, Mode of Action, and Antineoplastic Antibiotics

Communications of interest to microbiologists, chemists, pharmacologists, biochemists, clinicians, and pathologists

A valuable reference for laboratories and libraries

J. C. Sylvester, *Editor*

786 pages

\$15.00

1965

To: AMERICAN SOCIETY FOR MICROBIOLOGY
115 Huron View Blvd., Ann Arbor, Michigan

Please enter my order for — copy(s) of *Antimicrobial Agents and Chemotherapy—1964* at \$15.00 per copy. (Available in Europe from the H. K. Lewis Co., Ltd., London.)

I also wish to order:

_____ copy(s) of the 1963 volume at \$12.00 per copy.

_____ copy(s) of the 1962 volume at \$12.00 per copy.

_____ copy(s) of the 1961 volume at \$10.00 per copy.

☐ Bill me.


☐ Payment enclosed.

(Make checks payable to American Society for Microbiology.)

NAME _____

ADDRESS _____

new "slant" on an old friend

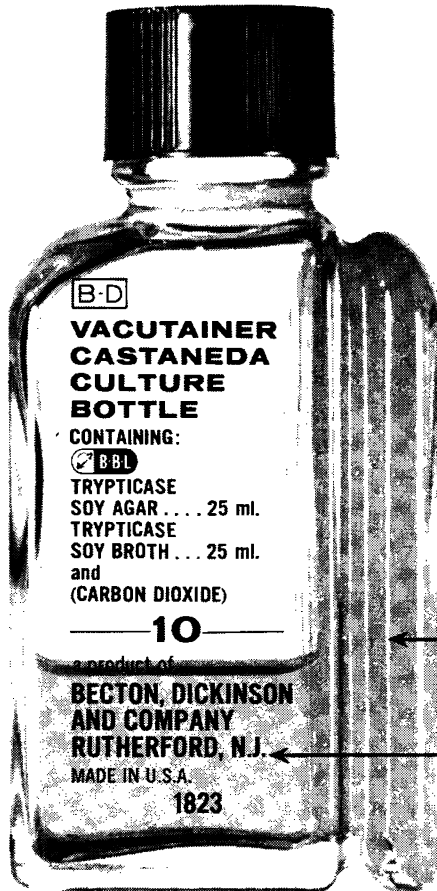

agar surface
and
fluid medium
in the
same bottle

B-D VACUTAINER CASTANEDA CULTURE BOTTLE

This new addition to the VACUTAINER Culture Bottle line provides both *surface* and *broth* growth—a 2-in-1 blood culture that yields more data, saves time and equipment.

As with all VACUTAINER Culture Bottles, contamination is avoided by blood collection through a closed, sterile system . . . additional bottles can be inoculated with the same venipuncture and B-D Blood Taking Unit . . . specimens can be taken *aerobically* or *anaerobically*.

a B-D  product

BECTON, DICKINSON AND COMPANY • RUTHERFORD, NEW JERSEY

In Canada: Becton, Dickinson & Co., Canada, Ltd., Toronto 10, Ontario

B-D, B-B-L, CASTANEDA, DISCARDIT, TRYPTICASE, AND VACUTAINER ARE TRADEMARKS.

10264



Focus on the finest



GOLD SEAL® SLIDES and COVER GLASSES

Synonymous with quality for almost half a century, microslides and cover glasses with the familiar GOLD SEAL label are guaranteed as perfect as the highest quality raw materials and skilled precision manufacturing can make them. *Continual research assures you of the best possible product for the best possible results in all applications.*

Gold Seal Slides come to you pre-cleaned, ready to use, in the *new* lint-free box. This box has an inner cellophane wrap as well as a sealed outer wrapper. The box also contains plastic inserts which readily convert the empty box into a Slide File convenient for transporting, filing and storing prepared slides. The special Stand-Rite tabs keep the slides upright in the box dur-

ing dispensing, permitting fingertip removal without smearing.

Gold Seal Cover Glasses are manufactured to the same exacting standards as Gold Seal Slides. They are made of a special, non-corrosive, non-fogging glass of uniform thinness. Each is individually selected and guaranteed perfect. Available in every convenient size and thinness, Gold Seal Cover Glasses are dispensed clean from lint-free plastic boxes.

Your dealer carries the Gold Seal Microslides and Cover Glasses in addition to a large selection of Clay-Adams microslide boxes, cabinets, and other useful accessories.

Clay-Adams
New York 10, N. Y.

SPORES III

A publication comprising the papers presented at the Third Spore Conference, held at Allerton Park, Ill., 2-4 October 1964.

Table of Contents

- W. G. MURRELL AND A. D. WARTH. Composition and Heat Resistance of Bacterial Spores.
VLADIMIR VINTER. Commencement of Synthetic Activities of Germinating Bacterial Spores and Changes in Vulnerability of Cells During Outgrowth.
MINORU UEHARA AND HILMER A. FRANK. Factors Affecting Alanine-Induced Germination of Clostridial Spores.
J. C. LEWIS, N. S. SNELL, AND GORDON ALDERTON. Dormancy and Activation of Bacterial Spores.
L. LEON CAMPBELL, CHRISTINA M. RICHARDS, AND ESTHER E. SNIFF. Isolation of Strains of *Bacillus stearothermophilus* with Altered Requirements for Spore Germination.
V. R. SRINIVASAN. Intracellular Regulation of Sporulation of Bacteria.
ROBERT W. BERNLOHR. Role of Amino Acids in Sporulation.
ARTURO FALASCHI, JAMES SPUDICH, AND ARTHUR KORNBERG. Deoxyribonucleic Acid Polymerase and Related Enzymes in Spores of *Bacillus subtilis*.
H. L. SADOFF, J. A. BACH, AND J. W. KOOLS. Significance of Multiple Forms of Glucose Dehydrogenase in Relation to Its Heat Resistance.
ROY H. DOI. Genetic Transcription of *Bacillus subtilis*.
JOHN SPIZZEN. Analysis of Asporogenic Mutants in *Bacillus subtilis* by Genetic Transformation.
I. TAKAHASHI. Mapping of Spore Markers on the Chromosome of *Bacillus subtilis*.
ERNST FREESE AND MICHAEL CASHEL. Initial Stages of Germination.
J. EDWARD DONNELLAN, JR., ELLA H. NAGS, AND HILLEL S. LEVINSON. Nucleic Acid Synthesis During Germination and Postgerminative Development.
JEKISIEL SZULMAJSTER AND RICHARD S. HANSON. Physiological Control of Sporulation in *Bacillus subtilis*.
A. KEYNAN AND HARLYN HALVORSON. Transformation of a Dormant Spore into a Vegetative Cell.
A. KEYNAN, G. ISSAHARY-BRAND, AND Z. EVENCHIK. Activation of Bacterial Spores.
J. C. VARY AND N. G. MCCORMICK. Kinetics of Germination of Aerobic *Bacillus* Spores.
Discussion. HILLEL S. LEVINSON AND MILDRED T. HYATT.
Y. KOBAYASHI, W. STEINBERG, A. HIGA, H. O. HALVORSON, AND C. LEVINTHAL. Sequential Synthesis of Macromolecules During Outgrowth of Bacterial Spores.
G. W. GOULD AND A. D. HITCHINS. Germination of Spores with Strange and Dark's Spore Lytic Enzyme.
HAROLD J. BLUMENTHAL. Changes in Carbohydrate Metabolism of *Bacillus cereus* Spores During the Return to Vegetative Growth.

L. Leon Campbell and H. Orin Halvorson, *Editors*

242 pages

Published 15 April 1965

\$6.00*

American Society for Microbiology
115 Huron View Blvd.
Ann Arbor, Mich. 48103

* Members of the American Society for Microbiology are entitled to purchase one copy at the member price of \$3.



Detection of the

SALMONELLA-SHIGELLA

Group

BACTO-S S AGAR

is a selective medium especially designed for use in isolation of *Shigella* and *Salmonella*. The selective action of this medium restrains the development of coliform bacteria and permits growth of the enteric pathogens. Because of the inhibitive action of the medium on coliform bacteria, it is possible to inoculate the medium heavily with feces thereby greatly increasing the number of positive isolations from samples containing very few pathogens.

BACTO-BISMUTH SULFITE AGAR

is a highly selective medium for isolation of *Salmonella*, especially *typhosa*. The unusual selective properties of this medium permit the use of large inocula of feces and other suspected material without overgrowth of extraneous intestinal bacteria.

BACTO-MacCONKEY AGAR

is an excellent differential medium for use in conjunction with Bacto-S S Agar and Bacto-Bismuth Sulfite Agar. This medium supports rapid and luxuriant growth of even the most fastidious strains of *Shigella* and *Salmonella*. Bacto-MacConkey Agar inhibits all Gram-positive organisms, permits growth of the coliform bacteria and affords excellent differentiation between the coliforms and lactose non-fermenting enteric organisms.

BACTO-BRILLIANT GREEN AGAR

is a highly selective medium recommended for the isolation of *Salmonella*. This medium permits rapid and luxuriant growth of *Salmonella* in 24 hours and at the same time inhibits other organisms encountered in stools, sewage, food products, tissue and other materials suspected of containing these microorganisms.

BACTO-TETRATHIONATE BROTH BASE and BACTO-SELENITE BROTH

are preliminary enrichment broths for the *Salmonella* group. These enrichment media are particularly valuable in detecting *Salmonella* in specimens containing small numbers of these organisms.

DIFCO LABORATORIES

DETROIT 1, MICHIGAN