



A) Skeletal Muscle Tissue Digest

FMO

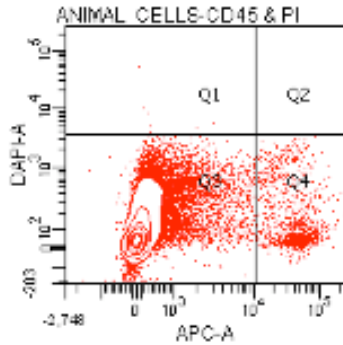

Tube: CD45 & PI			
Population	#Events	%Parent	%Total
■ All Events	50,000	####	100.0
■ P1	44,178	88.4	88.4
□ Q1	6	0.0	0.0
□ Q2	0	0.0	0.0
□ Q3	42,725	96.7	85.4
□ Q4	1,447	3.3	2.9


Stained for Endoglin

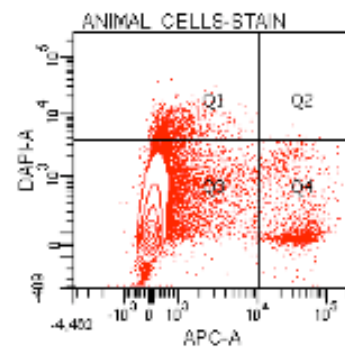

Tube: STAIN			
Population	#Events	%Parent	%Total
■ All Events	50,000	####	100.0
■ P1	43,469	86.9	86.9
□ Q1	818	1.9	1.8
□ Q2	8	0.0	0.0
□ Q3	41,310	95.0	82.6
□ Q4	1,333	3.1	2.7


B) Cultured Muscle Cells

FMO

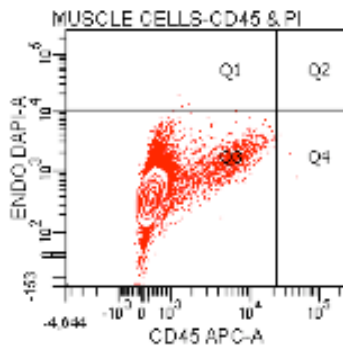

Tube: CD45 & PI			
Population	#Events	%Parent	%Total
■ All Events	50,000	####	100.0
■ P1	36,819	73.6	73.6
□ Q1	3	0.0	0.0
□ Q2	0	0.0	0.0
□ Q3	36,811	100.0	73.6
□ Q4	5	0.0	0.0

Stained for Endoglin

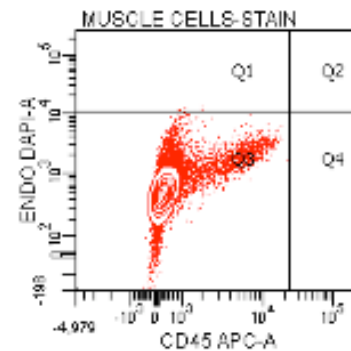

Tube: STAIN			
Population	#Events	%Parent	%Total
■ All Events	50,000	####	100.0
■ P1	36,608	73.2	73.2
□ Q1	6	0.0	0.0
□ Q2	0	0.0	0.0
□ Q3	36,601	100.0	73.2
□ Q4	1	0.0	0.0

Additional File 5: Specificity of endoglin antibody for endothelial cells in skeletal muscle digest compared to the lack of non-specific binding in a cultured muscle cell suspension.

While cell suspensions of collagenolytically digested skeletal muscle tissue from mice (A) show 1.9% of live cells staining for endoglin, 0% of C2C12 muscle cell suspensions (B) stain for endoglin demonstrating that non-specific binding of our antibody is not occurring. Fluorescence minus one (FMO) samples shown on the left lack endoglin antibody but are stained for CD45 and PI uptake. Endoglin+/CD45- cells are shown in Quadrant 1; Endoglin+/CD45+ cells are shown in Quadrant 2; CD45-/Endoglin- cells are shown in Quadrant 3; CD45+/Endoglin- cells are shown in Quadrant 4.