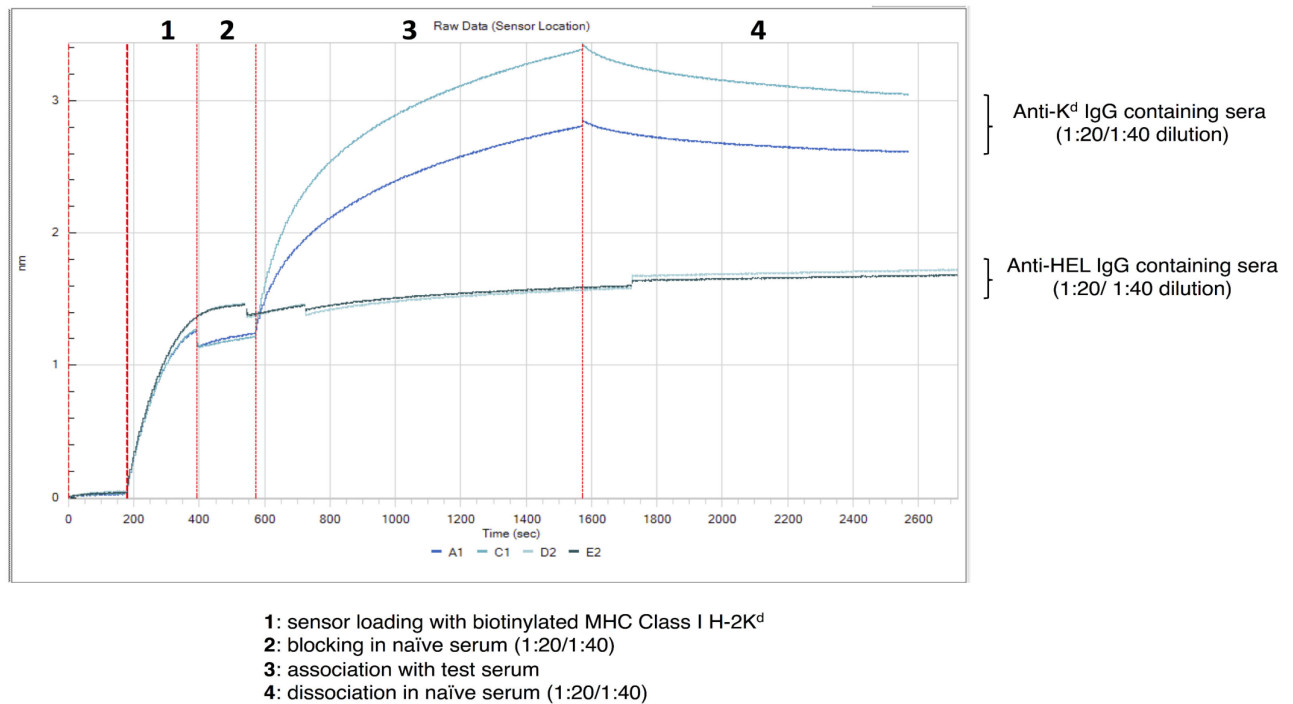


## SUPPLEMENTAL DATA

**Figure S1**



**Figure S1. Representative raw data binding kinetics of test sera containing H-2K<sup>d</sup> alloantibody.**

Streptavidin-coated sensors were moved between wells and binding kinetics were plotted in real time after baseline establishment and sensor pre-equilibration in the following order: **1)** loading with biotinylated H-2K<sup>d</sup> protein, **2)** 2<sup>nd</sup> baseline establishment and blocking with naïve BL/6 *Tcrbd*<sup>-/-</sup> serum (diluted at 1:20 or 1:40), **3)** association with test serum diluted at 1:20 or 1:40 (shown here from day 50 wild-type help-limited *Tcrbd*<sup>-/-</sup> recipients [top two curves] or from BL/6 SW<sub>HEL</sub> recipients challenged with hearts from BL/6.mHEL mice [bottom two curves]), and **4)** dissociation in naïve serum at a dilution equivalent to test serum. Sonograms were corrected for inter-step alignment on Y-axis, and baseline drift was corrected by subtracting the average value obtained from reference sample (H-2K<sup>d</sup> loaded sensor exposed to naïve B6 *Tcrbd*<sup>-/-</sup> serum) and reference sensor (unloaded sensor [no H-2K<sup>d</sup> protein] exposed to naïve B6 *Tcrbd*<sup>-/-</sup> serum).