


Supplementary Materials

Measurement of cellular ROS

NCI-H292 cells were plated into 96-well plates and grown overnight. Cells were washed with OPTI-MEM for two times and incubated with 3 μ M of CM-H2DCFDA (Invitrogen, Carlsbad, CA) for 30 min. Cells were then treated with 7.5 μ g/ml AFE in OPTI-MEM. The fluorescence was detected every 5 min by a plate reader (Tecan Infinite® M1000, Piedmont, NC). Background reading from the cells that were not loaded with CM-H2DCFDA was used as a blank.

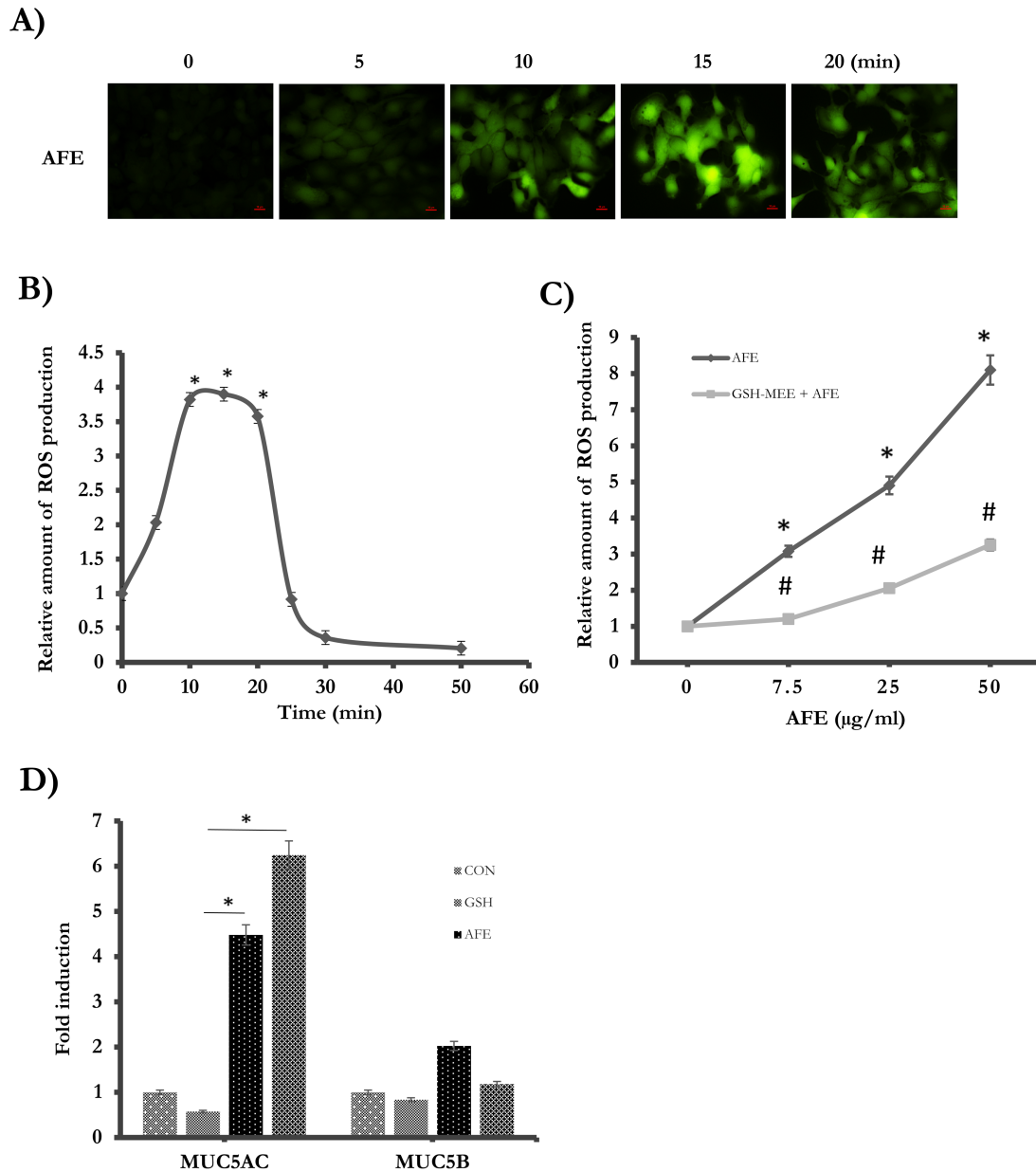

S2 Fig. (A) NCI-H292 cells were stimulated with 7.5 µg/ml AFE and intracellular ROS generation was measured every 5 min. (B) Quantification of ROS generation. $n=5$. (C) Application of GSH-MEE significantly blocked the dose-dependent AFE-induced ROS generation. (D) The cells were pre-treated with 5 mM GSH for 1 hr, and then treated with AFE for 6 hrs. MUC5AC and MUC5B were quantified by Real-Time PCR. *: AFE vs control (0 µg/ml) $P < 0.05$. #: AFE vs GSH-MEE + AFE, $P < 0.05$. $n=4$.